


BURN CAGE



3,500/-

Catheter Hook For Cardiac OT-Lab

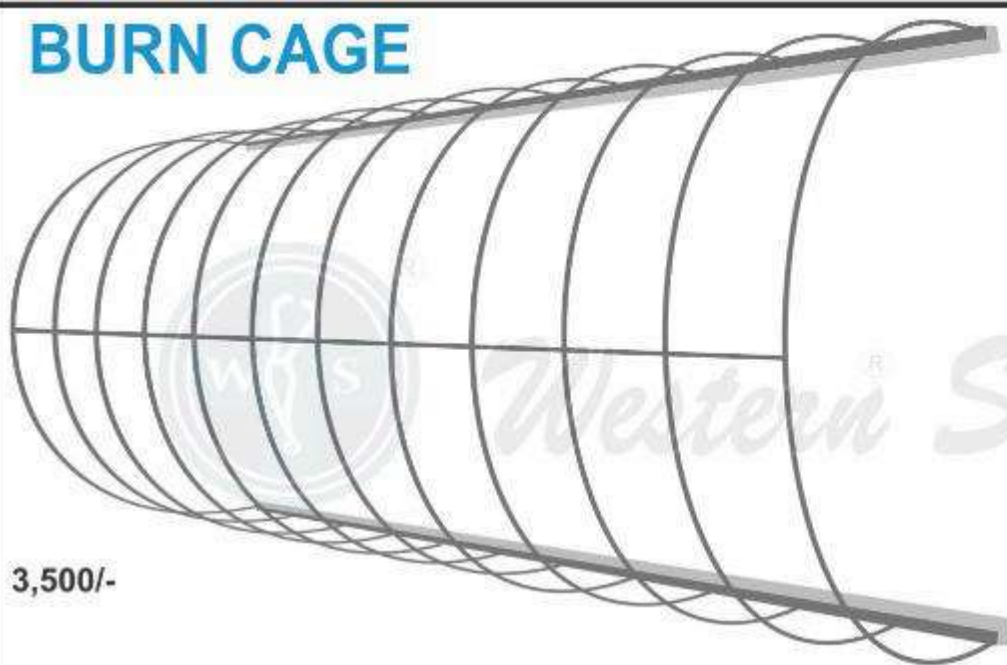

3,500/-

Variable Height Tilt Table




Height adjustable tilt table with two motors. One motor ensures effortless and smooth tilting; the other adjusts height and facilitates patient transfer. The Variable Height Tilt Table is the ideal standing aid for neuromotor revalidation and adaption exercises from lying to standing.

Another area where it can be applied is general orthostatic examination, more specifically examination and treatment of vasovagal syncope. Standard with three straps.

Standard model comes with a dual adjustable footplate. The tilt table also has a hoist access as standard and a quick release Trendelenburg positioning. It can also be fitted with an optional Torso Support for improved patient support and comfort, as well as an activity tray - for extra working space for hand activities.

Dimensions

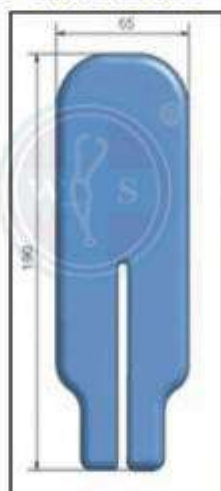

Height



Table Positions



Specifications

Minimum Height : 51cm/20"

Maximum Height : 86cm/33.9"

Cushion Width : 65cm/25.6"

Leg Length (+ Footplate) : 200cm/78.7"

Lift Mechanism : Electric

Breathing Hole : 0 / Wheels : 4 Braked

Tilting Mechanism : Electric

User Weight Rating : 205kg/451lb



* Carrying Warranty at Rajkot Office only on site warranty not provided free of cost.

* Due to constant upgradation design, price and features are subject to change any time without prior notice.



Western Surgical

WESTERN SURGICAL
EQUIPMENTS PRIVATE LIMITED

* www.westernsurgical.in *

Email : westernsurgical@gmail.com

Ph. +91-281-2232871 / 2243159

Mo. +91 94280 03750/52/54

© Copyright Reserved by : Western Surgical - Rajkot.